

# *ACE TURP LOOPS*

**Trans Urethral Resection of Prostate (TURP)** is a procedure done by urologists to reduce the size of the prostate gland, entering through the urinary tube (opening in the penis) of male patients. The prostate gland forms the opening at the end of male bladder, which allows the passage of urine from the bladder to the outside, through the Urethral canal in the penis. When the prostate gland becomes enlarged, passing of urine becomes difficult and it is then necessary to remove a portion of the gland so that the opening can be enlarged for proper urine passage. The enlargement can be done by cutting small slivers of the prostate or by vaporisation or both. Varying electrodes, inserted through the working element of a Resectoscope, are used for varying purposes.

**Cutting Loops :** is used to remove prostate tissue in the form of small chips.

**Collin's Knife :** is used to make straight incisions in a very hard prostate tissue.

**Roller Electrode :** is used to vaporize the prostate tissue using high power.

**Ball Electrode 3mm & 5mm :** is used to coagulate and smoothen the freshly resected surface of the prostate gland. 5mm electrode is used when the surface to be coagulated is large.

**Kriloktrode :** is a new electrode used for simultaneous cutting and vaporization.

**Cold Knife :** is used to make incisions of ureteral stenosis for upper urinary tract obstructions.

**Spike Roller 3mm & 5mm :** is used for electro - vaporisation in the treatment of BPH.




Design	Cutting	Colling's	Roller	Ball 3mm Ball 5mm	Kriloktrode (Patented)	Cold Knife	Necklace	Spike 3mm Spike 5mm
30° Storz 1 stem	W - SC	W - SZ	W - SR	W - SB W - SS	W - SK	W - SP	W - SN	W - SW W - SV
30° Storz 2 stem	W - DC	W - DZ	W - DR	W - DB W - DS	W - DK	W - DP	W - DN	W - DW W - DV
12° Olympus 1 stem	W - OC	W - OZ	W - OR	W - OB W - OS	W - OK			
30° Olympus 1 stem	W - LC	W - LZ	W - LR	W - LB W - LS	W - LK			
15° Wolf 1 stem	W - WC	W - WZ	W - WR	W - WB	W - WK			


- For 24 Fr loops add – 24F to Cat Nos. For 27 Fr Loops add – 27F to cat nos.



# Ace Nottingham Dilators

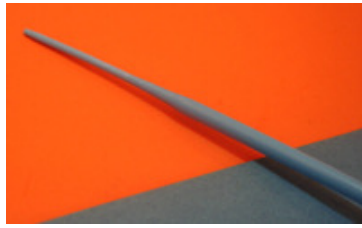

Used to dilate the intramural tunnel and ureter in one step. It may also be used to inject contrast medium to determine sub mucosal tunneling or anesthetic jelly to aid in dislodging ureteral stones. The 5cm distal taper allows gradual, single step dilation. Used with 0.038": guide wire. This device is Supplied sterile in a peel-open pouch. Intended for one-time use.

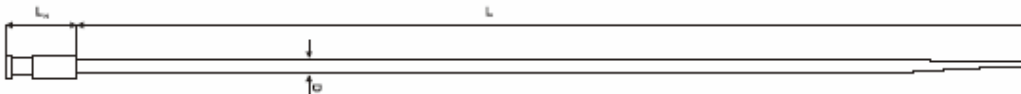

Sr No	Description	Cat Nos
01.	Nottingham Dilator 70 cm, 8/4 Fr	ND – 8/4
02.	Nottingham Dilator 70 cm, 9/5 Fr	ND – 9/5
03.	Nottingham Dilator 70 cm, 10/5 Fr	ND – 10/5
04.	Nottingham Dilator 70 cm, 10/6 Fr	ND – 10/6
05.	Nottingham Dilator 70 cm, 12/6 Fr	ND – 12/6
06.	Nottingham Dilator 70 cm, 14/6 Fr	ND – 14/6
07.	Nottingham Dilator 70 cm, 14/7 Fr	ND – 14/7
08.	Nottingham Dilator 70 cm, 16/6 Fr	ND – 16/6
09.	Nottingham Dilator 70 cm, 16/7 Fr	ND – 16/7
10.	Nottingham Dilator 70 cm, 16/8 Fr	ND – 16/8
11.	Nottingham Dilator 70 cm, 18/6 Fr	ND – 18/6

**\* Any other size can be manufactured on request.**

