




ESOPHAGEAL METALLIC STENTS

Used to maintain patency of malignant esophageal strictures and/ or to seal tracheoesophageal fistulas. Supplied sterile in a peel open pouch with deployer.

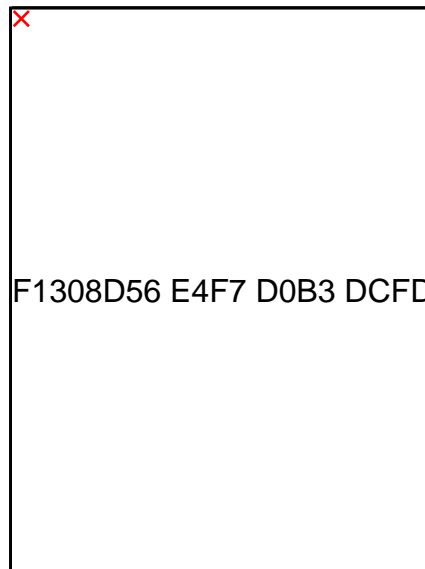

Key features:

- Available as Fully Covered, Both Ends Bare, With removal string
- Braided with single wire for optimum and uniform Radial Force and to maintain luminal patency
- High corrosion resistance and shape memory
- Slow expansion below body temperature
- Unique deploying mechanism is reconstrainable over 98% of deployment to aid in accurate repositioning of the stent
- Atraumatic stent endings to reduce tissue trauma
- Multiple radiopaque markers to facilitate accurate stent placement

Specifications

Size

Diameter 16,18,20,22,24,28 mm

Length 60, 80,100,120,140,150 mm