


HYBRID GASTRO GUIDEWIRE (STRIPED WITH 5.5 CM FLOPPY TIP)

Used to facilitate the placement of instruments during diagnostic or treatment procedures. Designed for single use only.

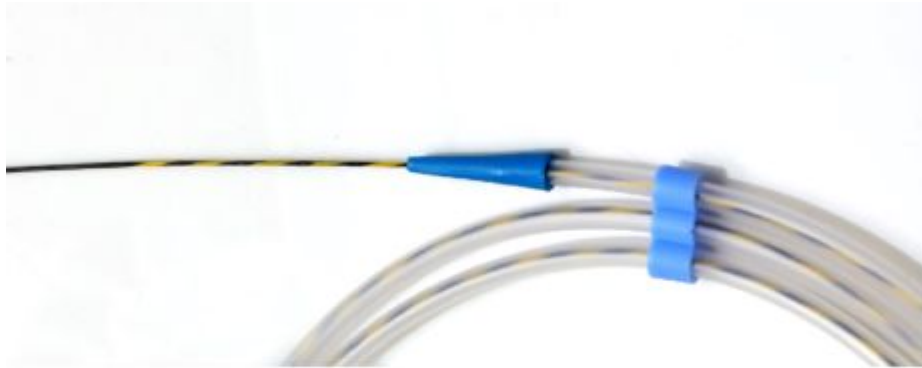

Key features:

- Straight, Flexible tip - for atraumatic advancement through gastrointestinal strictures
- Kink resistant Nitinol core with striped jacket - easy tracking of stents and instruments
- Low coefficient of friction for smooth advancement
- Better pushability for increased placement control
- 5.5 cm hydrophilic coated radiopaque tip provides smooth navigation with enhanced radiopacity
- Safe vessel access
- Available in Straight and Angled Tip

Specifications

Size

Diameter 0.018, 0.021, 0.025, 0.028, 0.032, 0.035

E - 1/7 Neptune Society, Kalyani Nagar, Pune - 411006, Maharashtra, India. Tel +91 20 26681483, Fax +91 20 26681500

Email info@acemedicaldevices.com acemedicaldevices.com

Inch

Length 260 and 450 cm