




## NOTTINGHAM ONE STEP DILATOR

Used to dilate the intramural tunnel and ureter in one step. May also be used to inject contrast medium to determine sub mucosal tunneling or anesthetic jelly to aid in dislodging ureteral stones. Designed for single use only.

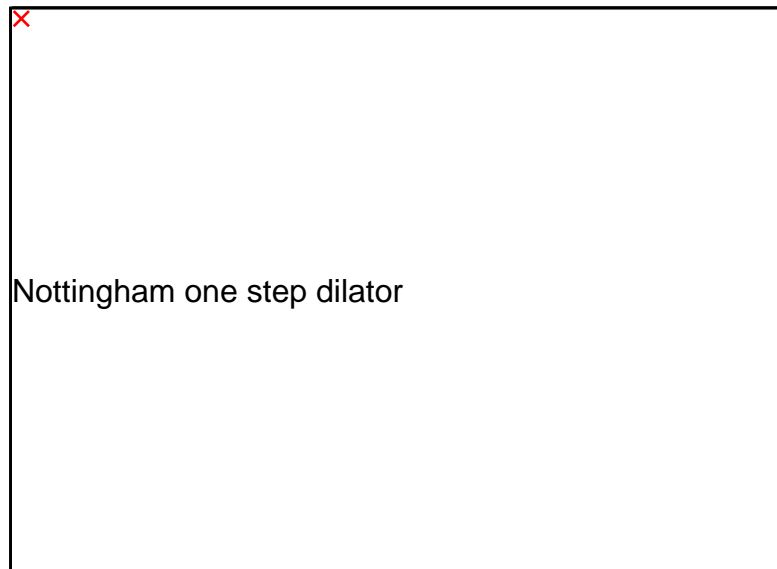

### Key features:

- Smooth surfaces for ease of introduction
- Uniform taper to reduce trauma.
- Used with 0.038" guide wire.

### Specifications

#### Size

**Diameter** 8/4, 9/5, 10/5, 10/6, 10/7, 12/6, 12/7, 14/6, 14/7, 16/6, 16/7, 16/8 and 18/6 Fr

**Length** 70 cm