


## PCN SET - PIGTAIL

PCN set is a complete set for percutaneous placement of a Pigtail Drainage Catheter in the renal pelvis for Nephrostomy drainage. Designed for single use only

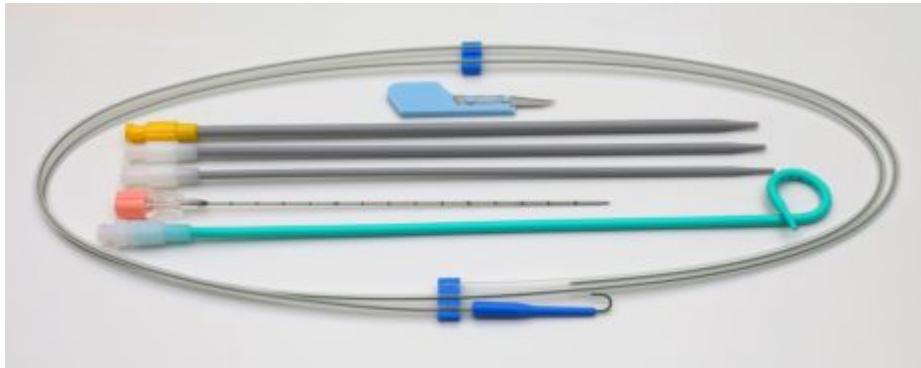

### Key features:

- Highly radiopaque for better visualization
- Pigtail tip configuration for the retention of the distal tip in the renal pelvis
- Looking type has more retention and can be kept for a longer duration
- Soft and flexible material for patient's comfort
- Available as PCN catheter with Lock
- Pigtail Catheter can be supplied with / without Hydrophilic coating

### Classification:

PCN Pigtail Set consists of:

- Pigtail Catheter
- Set of Fascial Dilators
- Initial Puncture Needle
- "J" Tip Guidewire
- Urine Bag Connector

E - 1/7 Neptune Society, Kalyani Nagar, Pune - 411006, Maharashtra, India. Tel +91 20 26681483, Fax +91 20 26681500

Email [info@acomedicaldevices.com](mailto:info@acomedicaldevices.com) [acomedicaldevices.com](http://acomedicaldevices.com)

- Scalpel

Available in Tray packaging as well as Pouch Packaging

## Specifications

### Size

**Diameter** 6, 7, 8, 8.5, 9, 10, 12, 14, 16,18  
Fr

**Length** 22 and 30 cm