


PCN DIRECT PUNCTURE SET- PIGTAIL

Designed for percutaneous Nephrostomy following the "Direct Puncture" technique or One Step Puncturing, particularly when there is a dilated urine collection system, a Hydro - Nephrosis.

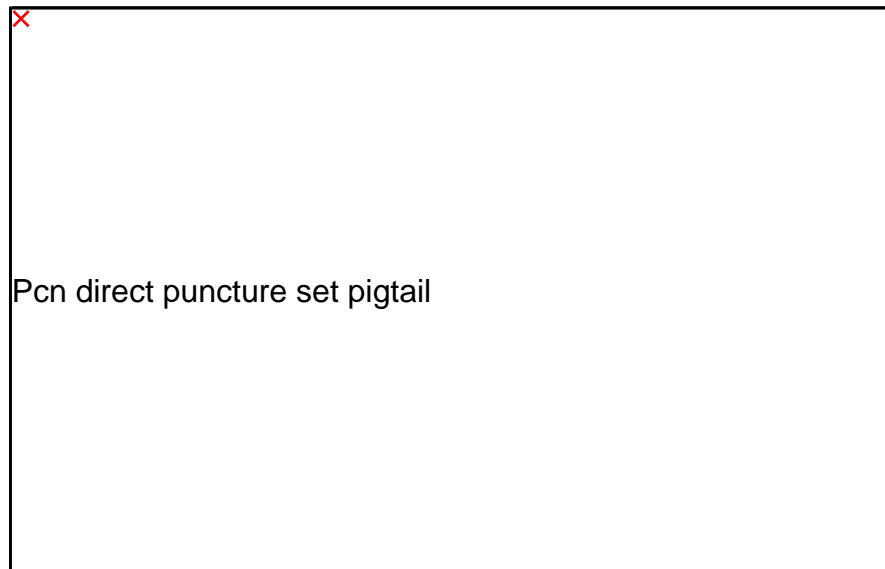

Key features:

- The catheter is back loaded on a two-part Trocar Needle and straightened and then introduced through a skin incision directly, using a Trocar Needle.
- Highly radiopaque for better visualization
- Pigtail tip configuration aids in the retention of the distal tip in the renal pelvis
- Locking type has more retention and can be kept for a longer duration.
- Soft and flexible material for patient's comfort

Classification:

Pigtail Catheter is supplied along with a Puncture Needle.

Available with / without Hydrophilic coating.

Specifications

Size

Diameter 6, 7, 8, 8.5, 9, 10, 12, 14, 16, 18 Fr

Length 22 and 30 cm