


## **RUNNER ENUCLEATION ELECTRODE - BIPOLAR**

Runner Electrode was developed as Transurethral enucleation with bipolar (TUEB) electrode. The electrode is supplied as a regular loop electrode along with a Runner, that easily helps isolate the prostate tissue from its capsule without much force. This allows for the en bloc enucleation of the whole adenoma in the prostate that was assumed to be difficult to treat. Product supplied Sterile and single use only.

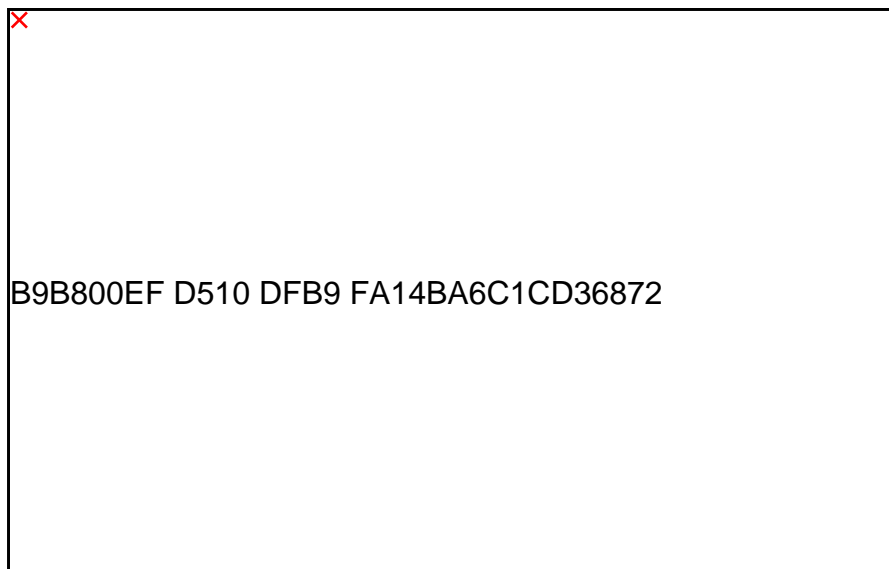

### **Key features:**

- Complete enucleation of the adenoma.
  - Better recuperation when compared to resection.
  - Minimum blood loss.
  - Better tissue sample for pathologic evaluation.

## Classification:

Compatible with

- STORZ 1 Stem Instrument
- Olympus 1 Stem Instrument
- Wolf 1 Stem Instrument

## Specifications

Size

**Diameter** 24 Fr